


Sphincterotomes

- Tip Orientation
Precise Positioning
- 2 & 3 Lumen Versions
Contrast Medium Injection while Guide Wire is in Place
- Monofilament Cutting Wire
Accurate Cut
- Version with Semi-insulated Cutting Wire
Safe & Controlled Papillary Roof Incision
- Color Marking with Black Rings
Macroscopically Controlled Cannulation Depth

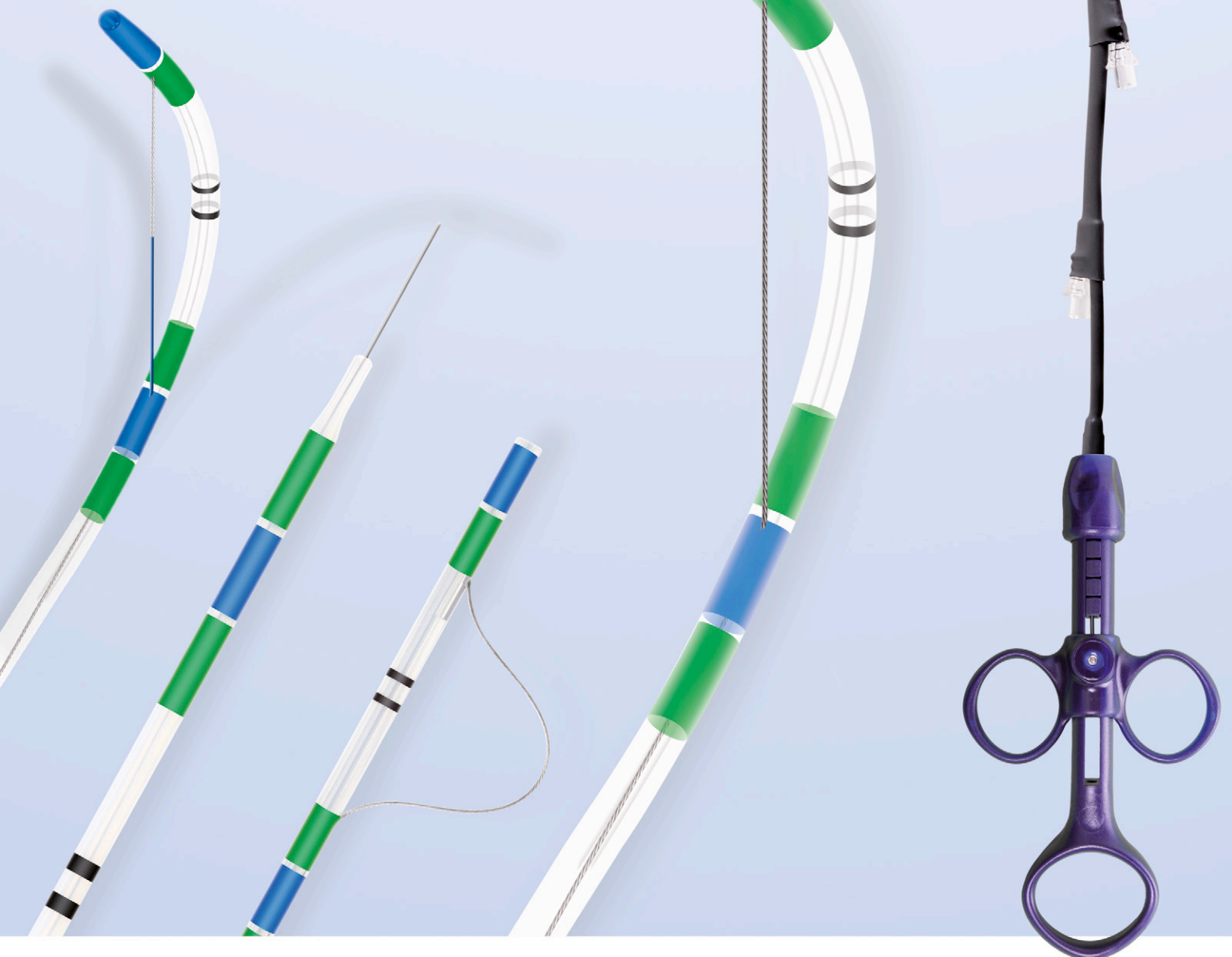

Sphincterotomes

Medi-Globe offers a wide range of sphincterotomes in triple, double and single lumen formats. Medi-Globe sphincterotomes have the ability to orientate the distal tip from 11 o'clock to 1 o'clock for precise positioning. Medi-Globe sphincterotomes are available with cutting wires of various lengths and in several tip styles to facilitate biliary and pancreatic duct access in nearly all anatomical requirements.

RotaCut®

RotaCut sphincterotomes are equipped with a braided cutting wire.

RotaCut® - M

RotaCut - M sphincterotomes are equipped with a monofilament cutting wire.

Premier Cut

Premier Cut sphincterotomes are equipped with a "semi-insulated", braided cutting wire.

S.O.S. Sphincterotomes

A S.O.S. "Speed Option Set" is comprised of a triple or double lumen sphincterotome, which is preloaded with a hydrophilic Teflon coated Nitinol guide wire in .035" diameter and a length of 450cm.



Medi-Globe sphincterotomes are directly compatible with Olympus electro-surgical generators. To verify compatibility with other electro-surgical generators please contact your Medi-Globe sales representative. Please also consult the Medi-Globe "Accessories" leaflet to see a complete list of active cords manufactured to fit most generators.

3-Lumen Guide Wire Sphincterotomes



Tapered Tip



RotaCut® – Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-31-20-020	Tapered Tip	5 mm	20 mm	2.55 mm	.035"
GSP-31-20-025	Tapered Tip	5 mm	25 mm	2.55 mm	.035"
GSP-31-20-030	Tapered Tip	5 mm	30 mm	2.55 mm	.035"

RotaCut® - M – Sphincterotome with Monofilament Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-32-20-020	Tapered Tip	5 mm	20 mm	2.55 mm	.035"
GSP-32-20-025	Tapered Tip	5 mm	25 mm	2.55 mm	.035"
GSP-32-20-030	Tapered Tip	5 mm	30 mm	2.55 mm	.035"

Premier Cut – Sphincterotome with Braided and Semi-insulated Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-27-20-030	Tapered Tip	5 mm	30 mm	2.55 mm	.035"

2-Lumen Guide Wire Sphincterotomes



Tapered Tip



RotaCut® – Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-21-20-020	Tapered Tip	5 mm	20 mm	2.55 mm	.035"
GSP-21-20-030	Tapered Tip	5 mm	30 mm	2.55 mm	.035"
GSP-24-20-020	Tapered Tip	10 mm (long nose)	20 mm	2.55 mm	.035"
GSP-24-20-030	Tapered Tip	10 mm (long nose)	30 mm	2.55 mm	.035"
GSP-21-17-020	Tapered Tip	5 mm	20 mm	1.75 mm	.021"
GSP-21-17-030	Tapered Tip	5 mm	30 mm	1.75 mm	.021"

2-Lumen Guide Wire Sphincterotomes



Precut



RotaCut® – Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-22-20-020	Precut Tip	-	20 mm	2.55 mm	.035"
GSP-22-20-030	Precut Tip	-	30 mm	2.55 mm	.035"

B II



RotaCut® – Billroth II Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-34-20-020	Billroth II	-	Variable	2.55 mm	.035"

Tapered Tip



Premier Cut – Sphincterotome with Braided and Semi-insulated Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Accepted Guide Wire
GSP-25-20-020	Tapered Tip	5 mm	20 mm	2.55 mm	.035"
GSP-25-20-030	Tapered Tip	5 mm	30 mm	2.55 mm	.035"

Pre-loaded Sphincterotome Sets

Tapered Tip



S.O.S. Sphincterotome – “Speed Option Set” with Pre-loaded Hydrophilically Teflon Coated Nitinol Guide Wire in .035” Diameter and 450 cm Length (GGW-03-35-450)

Catalog No.	Description	Tip Length	Cutting Wire	Outer Sheath Diameter	Accepted Guide Wire
GSP-32-20-031	3-Lumen RotaCut® – M – Sphincterotome, Tapered Tip	5 mm	Monofilament, 30 mm	2.55 mm	Pre-loaded
GSP-21-20-031	2-Lumen RotaCut® – Sphincterotome, Tapered Tip	5 mm	Braided, 30 mm	2.55 mm	Pre-loaded
GSP-25-20-031	2-Lumen Premier Cut – Sphincterotome, Tapered Tip	5 mm	Braided, Semi-insulated, 30 mm	2.55 mm	Pre-loaded

1-Lumen Sphincterotomes



Tapered Tip



Standard Tip



RotaCut® – Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter
GSP-01-17-020	Standard Tip	5 mm	20 mm	1.7 mm
GSP-01-17-030	Standard Tip	5 mm	30 mm	1.7 mm
GSP-04-17-020	Tapered Tip	5 mm	20 mm	1.7 mm
GSP-04-17-030	Tapered Tip	5 mm	30 mm	1.7 mm

Needle Knife



Needle Knife Sphincterotomes

Catalog No.	Needle Projection	Outer Sheath Diameter
GSP-11-17-005	0 - 5 mm	1.7 mm
GSP-11-17-010	0 - 10 mm	1.7 mm
GSP-11-17-015	0 - 15 mm	1.7 mm

Special Version

Fistulotome




Special Sphincterotomes

Catalog No.	Description	Outer Sheath Diameter
GSP-12-18-004	Fistulotom, Cutting Spiral 4 mm	1.85 mm

RotaCut®

Sphincterotomes – Single Use

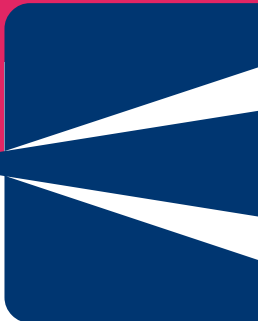

Distributed by:

L-GSP-SU-00-ENG Ver. 7.0 (2014-10)

Manufactured by: